



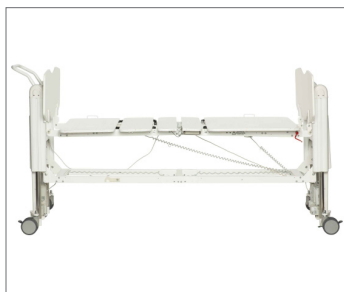

Multi Position  
Technology Ready  
Acute Care

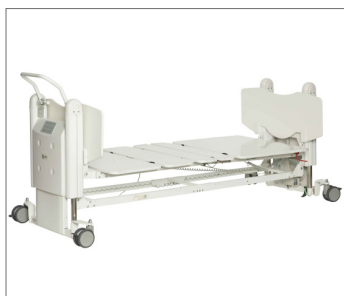

Our Floorline i Plus is a complete care management solution. With its extensive positioning, electronic CPR function and an open bus digital control box that allows for future upgrades, the Floorline i Plus floor level bed offers the highest level of care combined with the maximum versatility.

#### Clinical and Patient Benefits

The Floorline i Plus floor level bed is the perfect solution for preventing fall injuries in hospitals and acute care facilities. The bed lowers to a distance of 10cm/4" from the floor and raises to a working height of up to 80cm/31.5" to help promote existing safe patient initiatives and caregiver safety. At a distance of 10cm from the floor, the patient rolls out of the bed instead of falling, avoiding the dangerous and costly implications associated with a fall. In addition, the need for extra monitoring of these at risk patients is greatly decreased, adding to overall financial benefits for the hospital as less staffing is required.

#### Optimal Features

The Floorline i Plus features standard positions

such as Trendelenburg, reverse Trendelenburg and electronic cardiac positioning. With a full assortment of accessories, a built-in technology that allows for future upgrades, and durable LINAK electronics, the Floorline i Plus is a complete care management solution.

#### Quality

Our products are designed and manufactured to meet all the requirements of relevant standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone: +46 8 665 35 00  
Fax: +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street  
Port Melbourne, VIC 3207  
Phone: +61 1300 132 736  
Fax: + 61 (03) 9645 4954  
info.au@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com


#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories

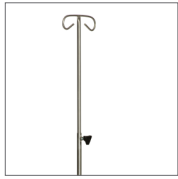

**IV Pole**  
Adjustable height. Can be mounted on either end of bed. Maximum SWL 7kg/15.4lb. Part no. IVP



**Care Assist Rail**  
Designed to assist patients safely in and out of a bed. Ergonomic design. Part no. 90702



**Extension Kit**  
Extends bed by 17.5cm/7" for taller patients. Part no. 90701



**Handset**  
Five function - included with bed. Part no. 75746  
Handset holder Part no. LC2250, Lock Out Key Part no. 75790



**Self Help Pole**  
Assists patient moving within confines of bed. Maximum SWL 75kg/165.3lb. Part no. 90700



**Side Rails**  
Available in 1/4 size [Part no. 72790], 1/2 size [Part no. 72791], Split 1/2+1/4 size [Part no. 72792], and 3/4 size [Part no. 72793]



**Wall Bumper Bar**  
Protects the head of the bed. Mounted on lower support bracket between wheels. Part no. 90706



**Oxygen Holder**  
Carries C size tank. Attaches to legs. Part no. 90707

## Features




**Push Handle**  
Push handle can be mounted on the foot end of the bed for easy operation. Replacement Part no. 90708



**Standard Positions**  
Includes auto-contour, backrest, kneebreak, Trendelenburg, reverse Trendelenburg, cardiac.



**Premium Casters**  
One central locking and three steering/directional casters for premium maneuverability.

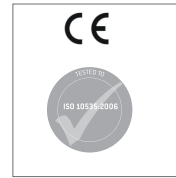

Manufactured to all relevant standards including IEC 60601-2-52:2009, EN1970:2002, AS/NZS 3200.2.38:2007

## Technical Specifications

Floorline i Plus			
Part no.	90500 (EU), 90501 (AU-1), 90502 (AU-2), 90503 (NA), 90504 (UK)		
Platform Position (Min-Max)	10cm - 80cm / 4" - 31.5"	Overall Bed Length (Standard to Extended with Extension Kit)	232.5cm - 250cm / 91.5" - 98.5"
Platform Dimensions (L x W)	200cm x 90.5cm / 78.75" x 35.5"	Platform Dimensions with Extension Kit (L x W)	217.5cm x 90.5cm / 85.5" x 35.5"
Sleeping Dimensions (L x W)	200cm x 90cm / 78.7" x 35.5"	Sleeping Dimensions with Extension Kit (L x W)	217.5cm x 90cm / 85.5" x 35.5"
Output Voltage	Max. 24V DC	IP Rating	IPX4
Input Voltage	120V, 60Hz (90404 to 90607), else 230V, 50Hz	Product Weight	150kg / 330lb
Platform Panel Angles	Backrest 70°, Thigh 45°, Knee 110°, Calf 25°, Trendelenburg/Reverse Trendelenburg 18°	Safe Working Load (SWL)	250kg / 550lb
Duty Cycle	2 minutes operation, 18 minutes rest	Maximum User Weight	185kg / 407lb

This product is manufactured in accordance with all relevant ISO standards.



**strack** Patienten- und Pflegebedarf

Strack AG 052 630 31 60  
Ebnatstrasse 125 www.strack.ch  
8200 Schaffhausen info@strack.ch

Learn More About This Product

### Follow Us on Social Media

