


strack Liftsysteme

FEATURES

- ◆ UL CERTIFIED ◆ 500 LB. (226 kg.) WEIGHT CAPACITY ◆ SETBACK RANGE 9" TO 36" (23 - 91 cm) ◆ CLEARANCES VARY BY WALL HEIGHT & WIDTH ◆ WATER DRAFT 10" ◆ 360° ROTATION ◆ FLIP UP ARMREST ◆ ALL STAINLESS STEEL CONSTRUCTION ◆ DURABLE POWDER COAT FINISH ◆ RECHARGEABLE BATTERY & WALL MOUNT CHARGER ◆ FIVE YEAR STRUCTURAL WARRANTY ◆ FIVE YEAR PRO-RATED ELECTRONICS WARRANTY

CLICK HERE TO VIEW THIS PRODUCT:
<https://aquacreekproducts.com/revolution/>

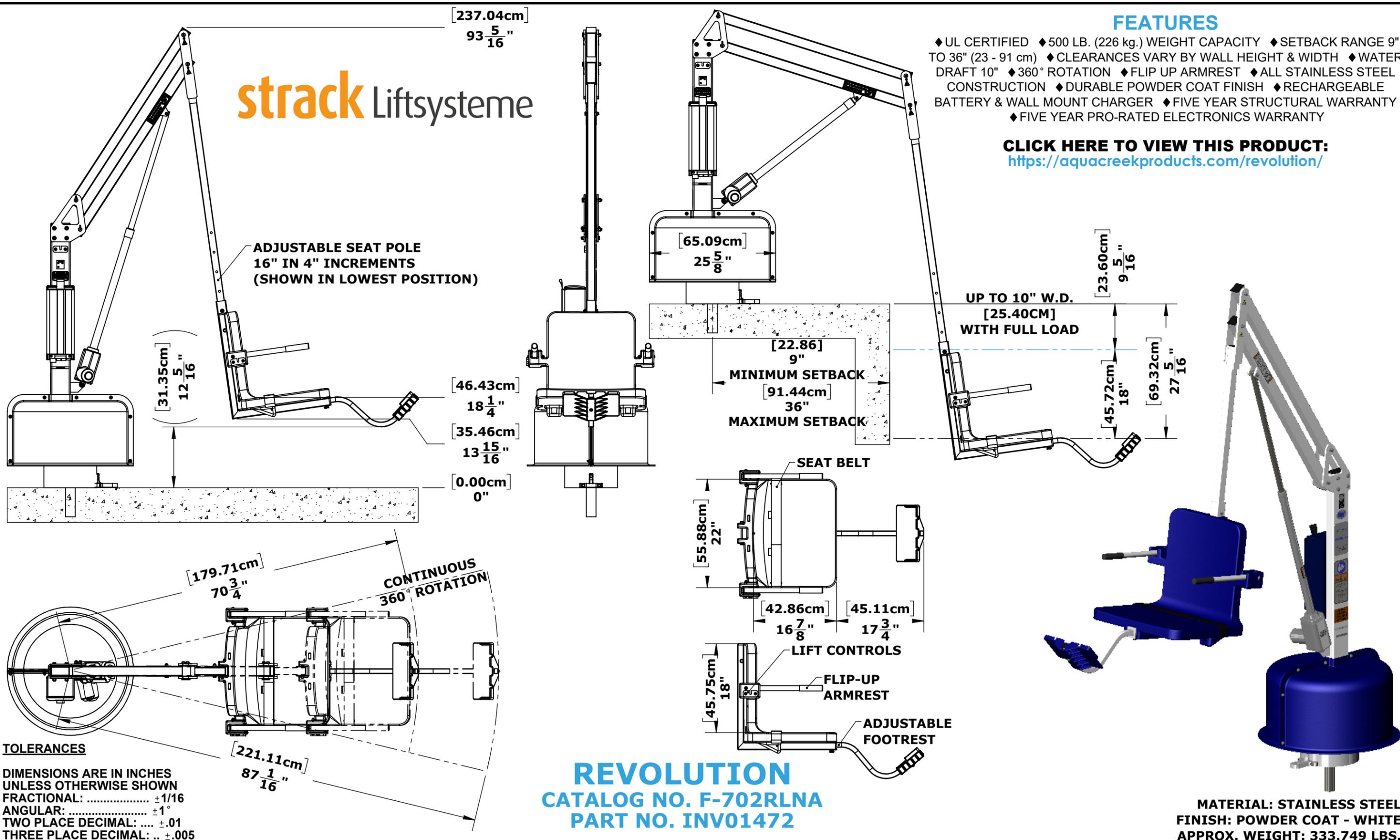

REVOLUTION - INV01472

REVOLUTION
CATALOG NO. F-702RLNA
PART NO. INV01472

MATERIAL: STAINLESS STEEL
FINISH: POWDER COAT - WHITE
APPROX. WEIGHT: 333.749 LBS.

TOLERANCES
DIMENSIONS ARE IN INCHES
UNLESS OTHERWISE SHOWN
FRACTIONAL: ±1/16
ANGULAR: ±1°
TWO PLACE DECIMAL: ±.01
THREE PLACE DECIMAL: .. ±.005

PROPRIETARY AND CONFIDENTIAL INFORMATION

THIS DRAWING AND ALL INFORMATION CONTAINED HEREIN IS THE SOLE PROPERTY OF AQUA CREEK PRODUCTS LLC. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT WRITTEN PERMISSION FROM AQUA CREEK IS PROHIBITED.

*SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT NOTICE

REVISIONS			
REV.	DESCRIPTION	DATE	APPROVED
P1	PRODUCTION RELEASE	4/19/2021	J.MAXWELL



9889 GARRYMORE LANE
MISSOULA, MT 59808
(406) 549-0769
FAX (406) 549-2602
www.AquaCreek.com

PRODUCT NAME: Revolution™ Pool/Spa Access Lift

CSI SECTION: 13 11 46 Swimming Pool Accessories

APPLICATIONS:

- In Ground Pools & Spas with up to **10-inch (25.40cm)** high by **22-inch (55.88cm)** wide curbs.
- On Ground & Partially Recessed Pools & Spas with walls up to **26-inches (66.04cm)** high by **4-inches (10.16cm)** wide when using the Optional pull out leg rest.

SUBMITTALS:

- Alternate make or model lifts will not be considered unless equal to the specified product in every respect. Alternative products must be submitted for approval no less than ten (10) calendar days prior to the bid due date. Submittal data must include complete documentation relating to all the specified features and include all manufacturer’s sales literature, specification sheets, installation/operation manuals, engineered drawings and documents that support the strength, service capacity, and the ability to meet the appropriate Federal Regulatory Standards associated with the product.

QUALITY ASSURANCE:

- The product must be a standard product of a company registered with the United States Food and Drug Administration (FDA) as a class 1 medical device manufacturer and that is regularly engaged in the manufacture of commercial grade swimming pool and spa accessibility equipment.
- The product shall be structurally capable of supporting a **500 pound (226.80kg)** load capacity and static load tested to 1.5 times the rated load.
- The product shall be certified to be compliant with UL test standard 60335-2-1000 Standard of Household and Similar Appliances: Particularly Requirements for Electrically Powered Pool Lifts.
- The product shall be tested and confirmed to comply with design standards outlined in the United States ADA Accessibility Standards for Accessible Designs.

SPECIFICATIONS:

- The products main structure shall be constructed of type 304 stainless steel; no dissimilar metallic components shall be allowed. The product shall be finished with a minimum 4-mil thick, UV resistant powder coat finish. All fastening hardware shall be of type 316 stainless steel. Product shall be supplied with a Linak Actuator suitable for medical use, a separate wall mount battery charger, a hand-held controller, and a 24-volt rechargeable battery.
- The chair of the product shall be of a rigid one-piece, polyethylene, molded design. The minimum seat width shall be **22-inches (55.88cm)** and the seat back shall be a minimum height of **17.5-inches (44.45cm)**. The seat shall feature dual flip-up arms and supplied with a 2-inch (5.08cm) wide adjustable safety belt with hook and loop fasteners.
- The seat height of the product will be adjustable from a load height of **18-inches (45.72cm)** to **34-inches (86.36cm)** in **4-inch (10.16cm)** increments to accommodate raised pool and spa walls.
- The product shall be capable of accommodating a deck to water distance (water draft) of up to **10-inches (25.40cm)** and capable of an anchor set back distance of up to **36-inches (91.44cm)**.
- An adjustable and removable footrest that moves with the seat shall be supplied with the product.

WARRANTY:

- Product shall have a minimum 5 year structural and 5 year pro-rated electronic warranty.

LIFT PART NUMBERS:

Item No	Order No	Description
F-702RLNA	INV01472	Revolution Lift w/Blue Seat (Specify White, Gray or Tan Powder Coat)
F-702RLNA-T	INV01782	Revolution Lift w/Tan Seat (Specify White, Gray or Tan Powder Coat)
F-702RLNA-G	INV01783	Revolution Lift w/Gray Seat (Specify White, Gray or Tan Powder Coat)
F-702RLNA-W	INV01784	Revolution Lift w/White Seat (Specify White, Gray or Tan Powder Coat)
F-702RLNA-C	INV01929	Revolution Lift (Choose Powder Coat and Seat Colors)

ANCHORING OPTIONS:

Item No	Order No	Description
F-710RLA	INV01480	Revolution Lift Anchor for Standard Concrete Decks
FC-77RLA-C	INV01562	Revolution Lift Semi-Recessed Anchor for 6” Thick Concrete

PRODUCT OPTIONS:

Item No	Order No	Description
F-0904CVR	INV01990	Revolution Lift Seat Cover in Blue
F-720SSC	INV00573	Revolution Lift Protective Lift Cover in Blue
F-720SSC-T	INV00575	Revolution Lift Protective Lift Cover in Tan
F-720SSC-G	INV00574	Revolution Lift Protective Lift Cover in Gray
F-045SCH	INV01692	Revolution Lift Solar Charger
F-422PLH-2	INV01455	Revolution Lift Adjustable Head Rest
F-423CS	INV01457	Revolution Lift Chest Strap
F-714RTC	INV01482	Revolution Lift Transport Cart
F-734RSA	INV01836	Revolution Lift Spine Board Assembly
F-LRB	INV01950	Revolution Lift Pull Out Leg Rest Blue
F-LRT	INV01952	Revolution Lift Pull Out Leg Rest Tan
F-LRG	INV01954	Revolution Lift Pull Out Leg Rest Gray
F-LRW	INV01953	Revolution Lift Pull Out Leg Rest White
F-019CA	INV01379	Revolution Lift Cycle Attachment

Aqua Creek Products LLC
 9889 Garrymore Lane
 Missoula, MT 59808
www.aquacreekproducts.com
 Tel: 888-687-3552 / 406-549-0769
 Fax: 406-541-8958



**9889 GARRYMORE LANE
 MISSOULA, MT 59808
 (406) 549-0769
 FAX (406) 549-2602
 www.AquaCreek.com**

PROPRIETARY AND CONFIDENTIAL INFORMATION

THIS DRAWING AND ALL INFORMATION CONTAINED HEREIN IS THE SOLE PROPERTY OF AQUA CREEK PRODUCTS LLC. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT WRITTEN PERMISSION FROM AQUA CREEK IS PROHIBITED.

***SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT NOTICE**

© 2021 AQUA CREEK PRODUCTS. ALL RIGHTS RESERVED.