


FEATURES

- ◆ UL CERTIFIED ◆ 300 LB. (136 kg.) WEIGHT CAPACITY ◆ WALL HEIGHT RANGE 21" TO 68" (54cm-172cm) ◆ PULL-OUT LEG REST ◆ 360° ROTATION ◆ FLIP UP ARMREST ◆ ALL STAINLESS STEEL CONSTRUCTION ◆ DURABLE POWDER COAT FINISH ◆ RECHARGEABLE BATTERY & WALL MOUNT CHARGER ◆ FIVE YEAR STRUCTURAL WARRANTY ◆ FIVE YEAR PRO-RATED ELECTRONICS WARRANTY

CLICK HERE TO VIEW THIS PRODUCT:
<https://aquacreekproducts.com/spa-lift-elite/>

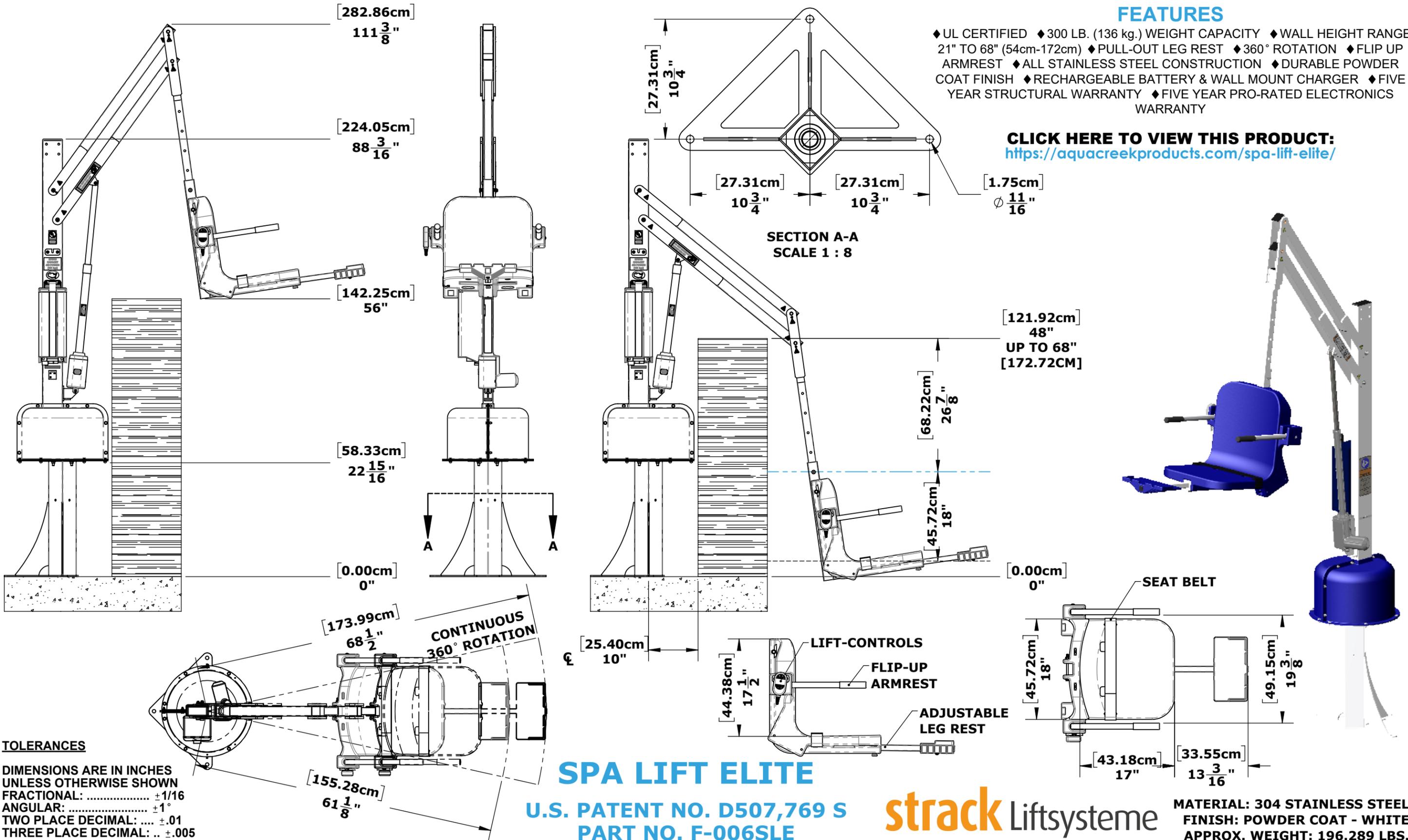

SPA LIFT ELITE

U.S. PATENT NO. D507,769 S
PART NO. F-006SLE

strack Liftsysteme

MATERIAL: 304 STAINLESS STEEL
FINISH: POWDER COAT - WHITE
APPROX. WEIGHT: 196.289 LBS.

TOLERANCES
DIMENSIONS ARE IN INCHES
UNLESS OTHERWISE SHOWN
FRACTIONAL: ±1/16
ANGULAR: ±1°
TWO PLACE DECIMAL: ±.01
THREE PLACE DECIMAL: .. ±.005

PROPRIETARY AND CONFIDENTIAL INFORMATION

THIS DRAWING AND ALL INFORMATION CONTAINED HEREIN IS THE SOLE PROPERTY OF AQUA CREEK PRODUCTS LLC. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT WRITTEN PERMISSION FROM AQUA CREEK IS PROHIBITED.

*SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT NOTICE

REVISIONS			
REV.	DESCRIPTION	DATE	APPROVED
P1	PRODUCTION RELEASE	4/29/2021	J.MAXWELL



9889 GARRYMORE LANE
MISSOULA, MT 59808
(406) 549-0769
FAX (406) 549-2602
www.AquaCreek.com

PRODUCT NAME: Spa Lift Elite Pool/Spa Access Lift

CSI SECTION: 13 11 46 Swimming Pool Accessories

APPLICATIONS:

- On Ground & Partially Recessed Pools & Spas with walls from **21-inches (53.34cm)** up to **68-inches (172.72cm)** high and up to **14-inches (35.56cm)** wide.

SUBMITTALS:

- Alternate make or model lifts will not be considered unless equal to the specified product in every respect. Alternative products must be submitted for approval no less than ten (10) calendar days prior to the bid due date. Submittal data must include complete documentation relating to all the specified features and include all manufacturer's sales literature, specification sheets, installation/operation manuals, engineered drawings and documents that support the strength, service capacity, and the ability to meet the appropriate Federal Regulatory Standards associated with the product.

QUALITY ASSURANCE:

- The product must be a standard product of a company registered with the United States Food and Drug Administration (FDA) as a class 1 medical device manufacturer and that is regularly engaged in the manufacture of commercial grade swimming pool and spa accessibility equipment.
- The product shall be structurally capable of supporting a **300 pound (136.08kg)** load capacity and static load tested to 1.5 times the rated load.
- The product shall be certified to be compliant with UL test standard 60335-2-1000 Standard of Household and Similar Appliances: Particularly Requirements for Electrically Powered Pool Lifts.
- The product shall be tested and confirmed to comply with design standards outlined in the United States ADA Accessibility Standards for Accessible Designs.

SPECIFICATIONS:

- The products main structure shall be constructed of type 304 stainless steel; no dissimilar metallic components shall be allowed. The product shall be finished with a minimum 4-mil thick, UV resistant powder coat finish. All fastening hardware shall be of type 316 stainless steel. Product shall be supplied with a Linak Actuator suitable for medical use, a separate wall mount battery charger, a hand-held controller, and a 24-volt rechargeable battery.
- The chair of the product shall be of a rigid one-piece, polyethylene, molded design. The minimum seat width shall be **18-inches (45.72cm)** and the seat back shall be a minimum height of **17.5-inches (44.45cm)**. The seat shall feature dual flip-up arms and supplied with a 2-inch (5.08cm) wide adjustable safety belt with hook and loop fasteners.
- The product shall be capable of accommodating an anchor set back distance of **10-inches (25.40cm)**.
- An adjustable and removable footrest that moves with the seat shall be supplied with the product.

WARRANTY:

- Product shall have a minimum 5 year structural and 5 year pro-rated electronic warranty.

LIFT PART NUMBERS:

Item No	Order No	Description
F-006SLE	INV01376	Elite Spa Lift w/Blue Seat (Specify White, Gray or Tan Powder Coat)
F-006SLE-T	INV01797	Elite Spa Lift w/Tan Seat (Specify White, Gray or Tan Powder Coat)
F-006SLE-G	INV01798	Elite Spa Lift w/Gray Seat (Specify White, Gray or Tan Powder Coat)
F-006SLE-W	INV01799	Elite Spa Lift w/White Seat (Specify White, Gray or Tan Powder Coat)
F-006SLE-C	INV01914	Elite Spa Lift (Choose Powder Coat and Seat Colors)

ANCHORING OPTIONS:

Item No	Order No	Description
F-SPAACK8	INV01543	Elite Spa Lift Anchor for Standard Concrete Decks
F-107UEWDA	INV01423	Elite Spa Lift Anchor for Wood Deck Applications

PRODUCT OPTIONS:

Item No	Order No	Description
F-0904CVR	INV01990	Elite Spa Lift Seat Cover in Blue
F-450SLECS	INV00544	Elite Spa Lift Protective Lift Cover in Blue
F-450SLECS-T	INV00546	Elite Spa Lift Protective Lift Cover in Tan
F-450SLECS-G	INV00545	Elite Spa Lift Protective Lift Cover in Gray
F-045SCH	INV01692	Elite Spa Lift Solar Charger
F-422PLH-2	INV01455	Elite Spa Lift Adjustable Head Rest
F-423CS	INV01457	Elite Spa Lift Chest Strap
F-019CA	INV01379	Elite Spa Lift Cycle Attachment

Aqua Creek Products LLC
9889 Garrymore Lane
Missoula, MT 59808
www.aquacreekproducts.com
Tel: 888-687-3552 / 406-549-0769
Fax: 406-541-8958



9889 GARRYMORE LANE
MISSOULA, MT 59808
(406) 549-0769
FAX (406) 549-2602
www.AquaCreek.com

PROPRIETARY AND CONFIDENTIAL INFORMATION

THIS DRAWING AND ALL INFORMATION CONTAINED HEREIN IS THE SOLE PROPERTY OF AQUA CREEK PRODUCTS LLC. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT WRITTEN PERMISSION FROM AQUA CREEK IS PROHIBITED.

*SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT NOTICE

© 2021 AQUA CREEK PRODUCTS. ALL RIGHTS RESERVED.

MODEL
P1
REVISION

SPA LIFT ELITE - F-006SLE

51-163-220

SCALE
1:18

SHEET
2