
EN
DE
FR

Instructions for Use
MOBIO Comfort Premium

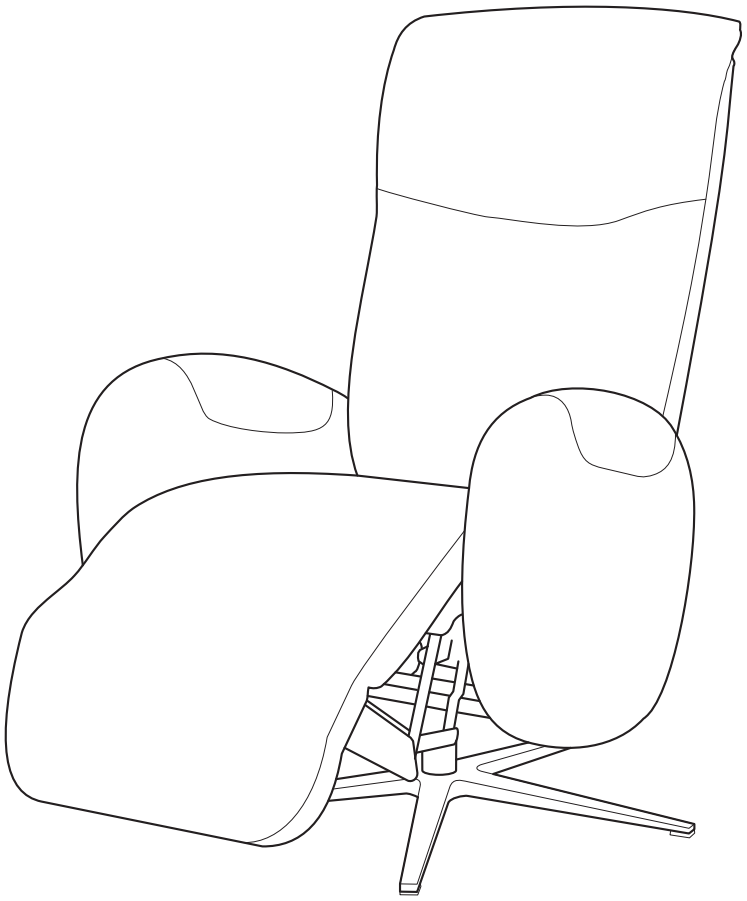2 YEAR WARRANTY

EN

Illustrations
Dimensions
Technical data
Article numbers

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

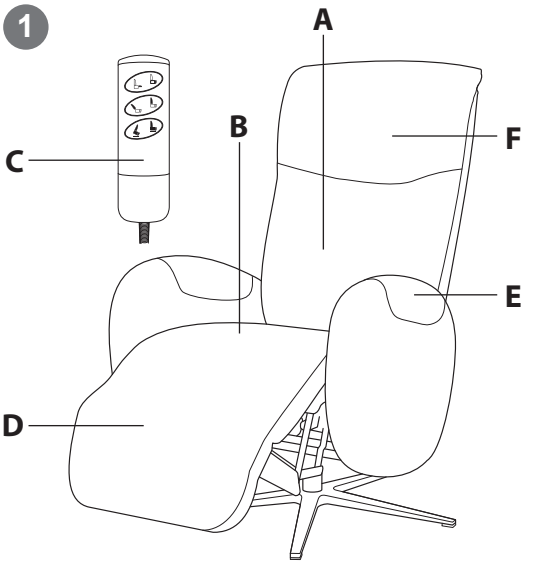
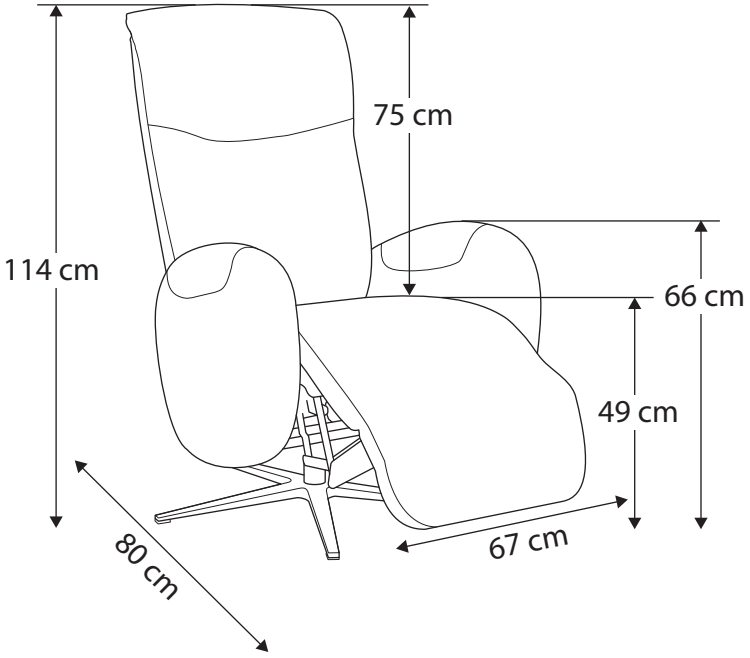
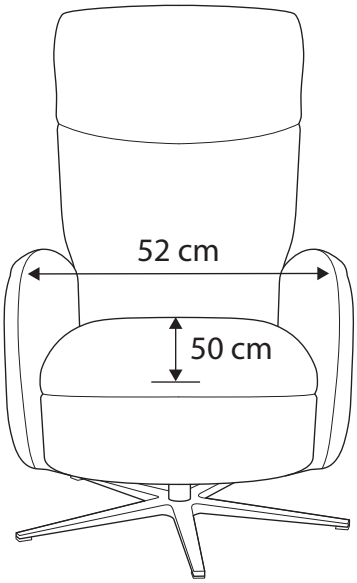
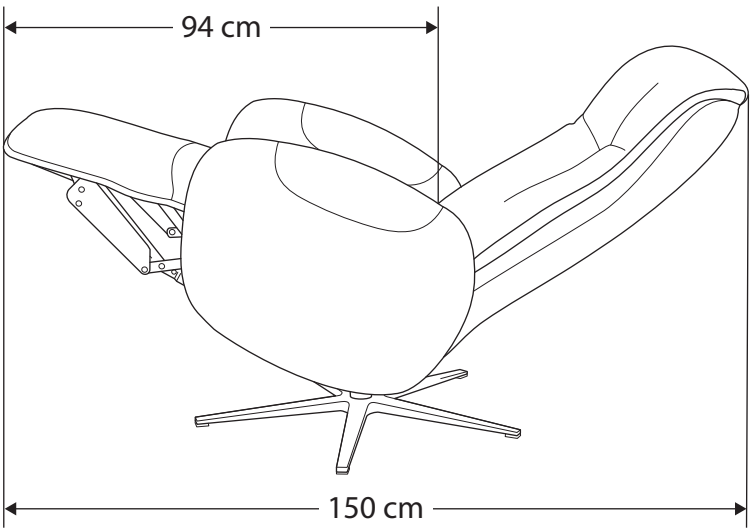
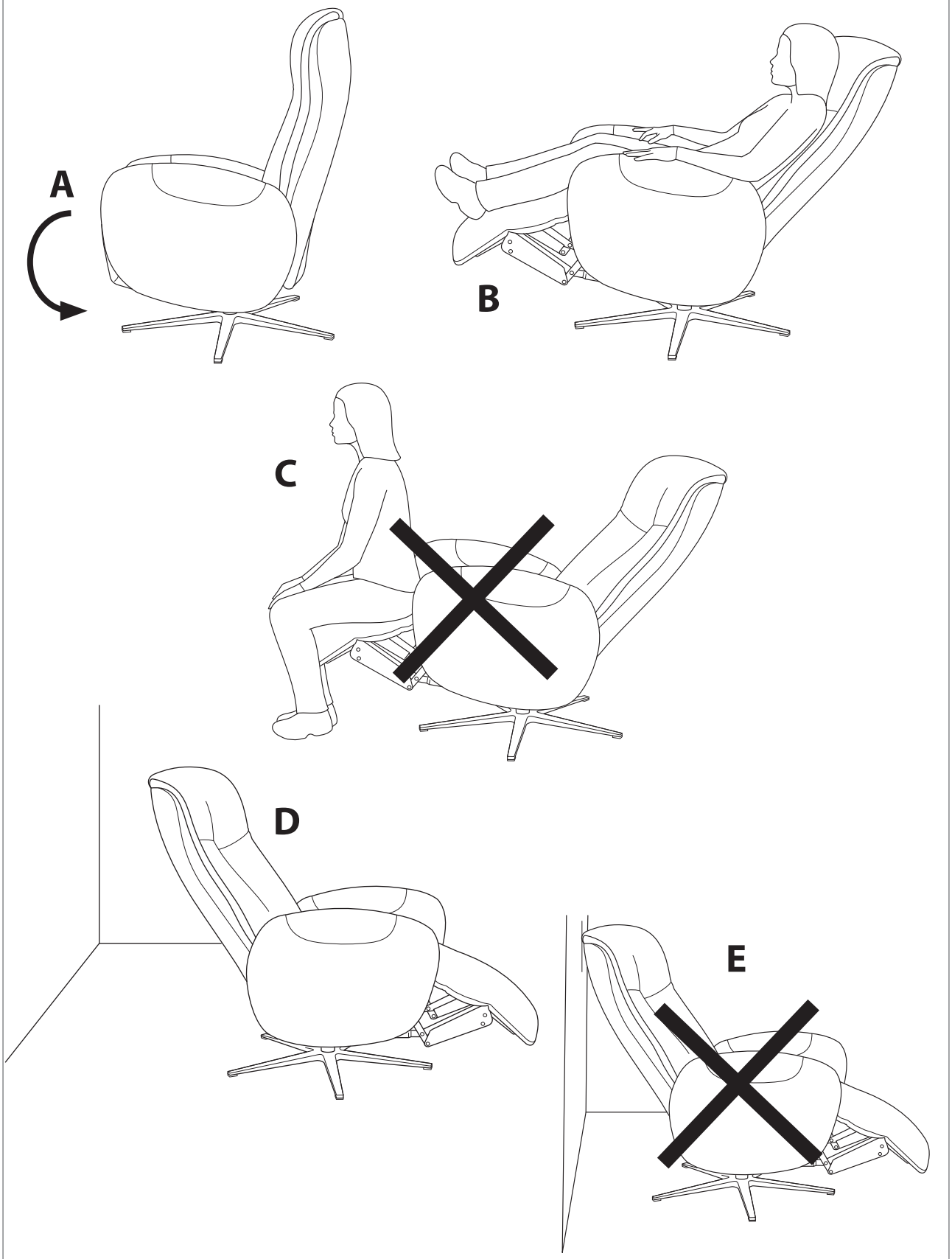
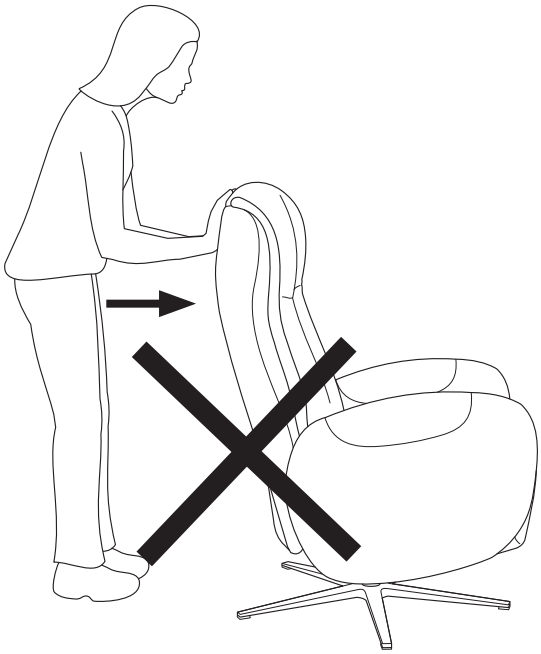
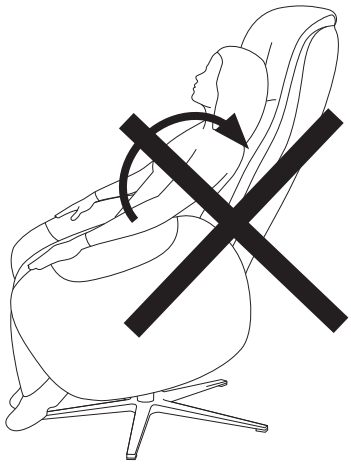
3

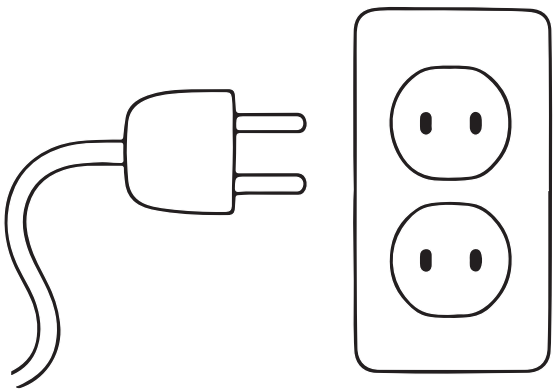4

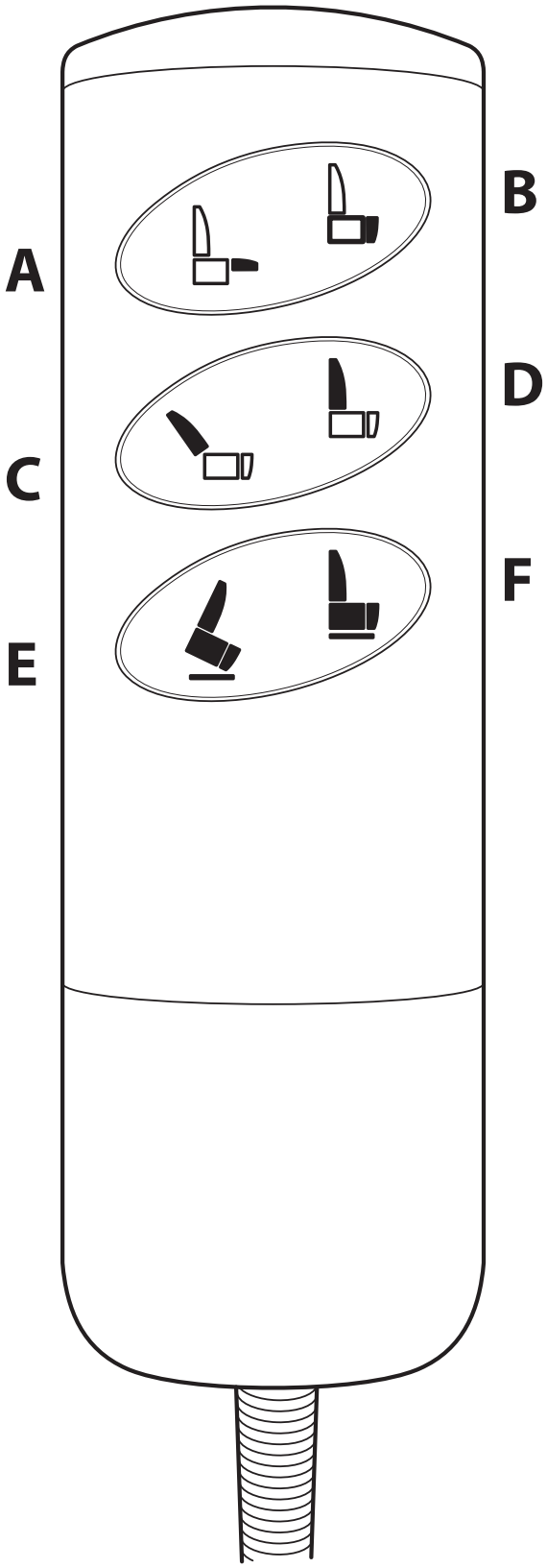Congratulations on your new MOBIO rise and recline chair

Our stand-up and reclining chairs have been developed with a special focus on your comfort. Thanks to ergonomically designed seats and adjustable footrests and backrests, they offer optimum support for your back and legs. MOBIO stand-up and reclining chairs are available in a variety of modern colours and attractive upholstery. They have sophisticated functions that increase comfort and safety in everyday life. Please check your product as soon as you receive it. If you have any questions, please contact your specialist mobility shop or dealer. If you have any difficulties reading this user manual, an electronic version is available at **mobio.be**.

Intended use

MOBIO Comfort Premium is intended for able-bodied adults, and the chair should not be equated with wheelchairs, medical chairs or beds for disabled people. The chair is intended for indoor use and can only be used by one person at a time within the weight limit of 160 kg. The user must be able to operate the product without external assistance. The warranty is void if the product is damaged by a user who does not comply with these restrictions.

0 The grey numbered circles in this user manual refers to the numbered information graphic and illustrations inside the cover

1 Get to know your MOBIO product

A Back rest
B Seat
C Remote control
D Footrest
E Armrests
F Head rest

2 Safety ⚠

Never enter or exit the chair without first returning the footrest to its original position (A). The footrest is designed only for its intended use (B) and should never be used as a seat or for similar purposes (C), as the hardware of the product may be damaged.

The chair should be positioned in such a way that no obstacles can impede the operation of the lifting and reclining mechanism (D). Under no circumstances should the backrest or footrest be placed against a wall (E), due to the risk of significant damage that is not covered by the warranty.

The chair must be placed on a level floor.

Do not push the chair forward (F).

Do not pull the chair sideways by its armrests (G).

Always take the precaution by disconnecting the power supply to the chair when cleaning or moving it, or when children are present unsupervised.

The chair is not a toy and should not be used by children. Be sure to use the chair with caution in the presence of children or pets.

Never recline the backrest (H) while the chair is in a raised position (I) due to the risk of injuries.

Do not pull the remote-control cord more than necessary. Do not sit on the remote control or let it fall to the ground. Do not allow children to use the remote control.

Carefully ensure that the product is properly installed and connected to the power supply before use.

Do not modify the product, as this may compromise your safety and void the warranty.

Check that the connectors under the chair are securely in place and not disconnected due to household cleaning, such as with a broom or vacuum cleaner.

Keep the chair away from light or heat sources.

Art. No. User Manual: 105145 – Revision A, 2025-09

3 Installation

Wired Power Supply
To activate the motors, plug it into the power outlet. Once connected to the socket, you can operate the remote control.

4 Adjustments MOBIO Comfort Premium

- **A** Rise the footrest
- **B** Lower the footrest + Rise support
- **C** Recline the backrest
- **D** Rise the backrest
- **E** Stand up support
- **F** Sit down support + Rise the footrest

ⓘ The neck support can be adjusted manually as required.

Troubleshooting

Problem	Probable Cause	Solution
The chair is not working	No voltage	Check that the battery is connected to the motor(s) and that the electrical outlet is working
	Defective chair	Contact customer service
The chair suddenly becomes unresponsive, no movement takes place	The heat protection has been triggered	Leave the chair for 20 to 30 minutes and try again
	A cable has been disconnected	Check that the remote-control cord is connected
	The battery fuse may be out of order	Contact customer service

Please contact your supplier for further assistance if everything above has been checked and the chair still does not function properly.

Storage

For long-term storage, measures must be taken to protect and preserve the chair:
• Protect the chair from sunlight to counteract the natural process of colour fading by covering the chair.
• Store the chair in dry and well-ventilated rooms.

Cleaning

Protect the chair from all sources of direct light (sun, moon, halogen lamps, etc.) to avoid the risk of discolouration.

Use a vacuum cleaner and regularly wipe the chair without excessive friction to avoid damaging the upholstery. The vacuum cleaner will remove dry dust but not greasy dust that may adhere to the upholstery. Therefore, it is recommended to clean the chair every six months. Clean the parts in direct contact with the scalp and hands with water and a neutral soap. Avoid household detergents or products containing solvents or oils that can damage the upholsteries. Always test on a hidden part of the chair.

Oily materials	Use a dry cloth to absorb the oil and remove the stain. If this is not sufficient, apply a small amount of dish soap and gently rub with a damp cloth.
Coffee, tea and beer	Clean with hot water. If the stain persists, dab it with a soapy solution.
Chewing-gum	Place an ice cube wrapped in cloth to cool the gum. Then scrape it off very gently.
Ballpoint pen ink	Dampen a cotton swab with a little vinegar and run it over the stain.
Food stain and blood	Carefully remove the superficial part and moisten the stain with some water. Then dab with a cloth soaked in soapy water.

Tip: For tough stains, consult a stain removal specialist.

Inspection / Maintenance

The mechanism should be dedusted regularly and the fastenings checked and retightened as required (at least once a year).The mechanism may exhibit slight play during use and minor noises or creaking noises cannot be ruled out during operation. In such cases, it is recommended to apply a small amount of lubricant to the joints. some surface coatings – especially artificial leather – can squeak during operation. The use of talcum powder is recommended to minimise these effects.

Material

	MOBIO Comfort Premium	
Seat	Metal frame with metal springs High resilience foam 34 kg/m³	Memory foam seat mattress 52 kg/m3, thickness 2.5 cm
Back rest	Metal structure High resilience foam 26 kg/m³ Manual articulated headrest	Memory foam seat mattress 52 kg/m3, thickness 2.5 cm
Footrest	Panel structure High resilience foam 34 kg/m³	Memory foam seat mattress 52 kg/m3, thickness 2.5 cm
Armrest	Plywood structure Polyether foam 20 kg/m³	
Swivel Multi-position lifting mechanism, 3 motors, ZERO-GRAVITY Stawett (made in germany)		

Warranty

The MOBIO Comfort Premium is guaranteed against defects in materials and workmanship for a period of 2 years, please contact your local specialist mobility shop or dealer in the event of a warranty claim, the warranty only applies if the wheelchair is used correctly in accordance with the instructions in this user manual.

Exclusions from Warranty

Unintentional damages caused during delivery or installation of the product by the distributor or user.

Damages caused by misuse of the chair.

Changes in appearance caused by normal use and wear of the product.

Products in abnormally dirty condition and/or with unintentional or deliberate damage.

Damages resulting from prohibited use: sitting on the footrest or armrests, entering or exiting the chair without first folding the footrest, moving the product by pulling sideways on the armrest or the backrest, obstructing the free movement of the elements in the backrest or footrest, especially against a wall or furniture.

Presence of a mark on any part of the chair indicating it was placed against a wall, especially on the backrest, will invalidate the product warranty.

For the electrical part, damages caused by power surges or lightning, poor electrical installation, cutting or tearing of electrical wires, breaking connection plugs, unintentional disconnection of motors or remote controls.

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MOBIO Switzerland GmbH